

# FLECOTEC

## VACUUM BAG

### VACUUM BAG SYSTEM




### VACUUM BAG for tablet presses

The special design of the vacuum bag has tabs attached on the outside, which hold the bag open in a frame. The flexible Flecozip MTS secondary interface in turn is located around the MTS transfer tube.

This arrangement ensures a closed entry of tablets into the flexible vacuum bag even under negative differential pressure conditions. The filled bag can be easily removed from the frame.

Whether standardized sizes or custom-made, the Vacuum Bag application provides the right connection and safe, contained transport from the tablet press to further processing units.

### COMPONENTS

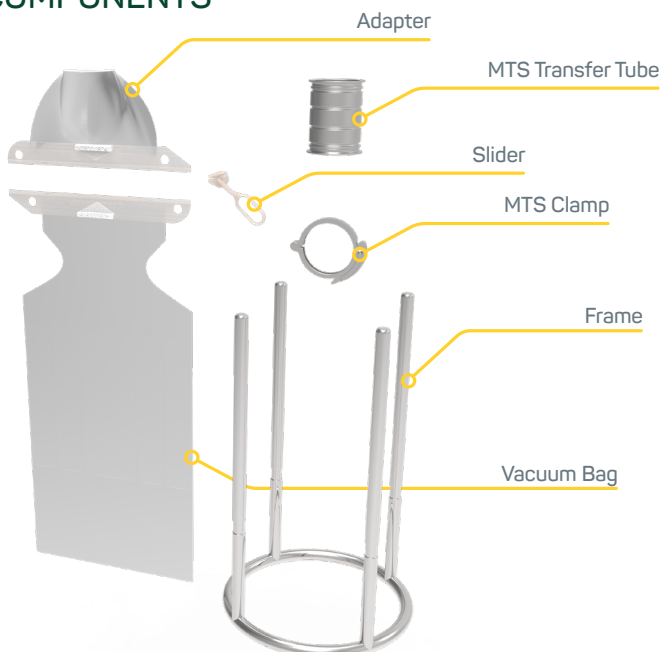



Flecotec AG

Breisgauring 13  
79427 Eschbach  
Germany

T + 49 7632 828 553-0  
F + 49 7632 828 553-9

[mail.ft@rommelag.com](mailto:mail.ft@rommelag.com)  
[www.rommelag.com](http://www.rommelag.com)

